


# Application 283-00

## Agilent Refinery Gas Analyzer

### Hydrocarbons in Extended Refinery Gas and Cracked Gas Analysis

### Technical Overview

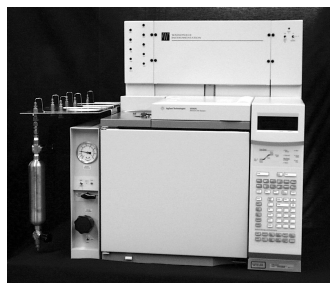

### Application Highlights


- A Flame Ionization Detector (FID) to detect the C1 through C7 paraffins and olefins to a lower detection limit of 20 ppm, except for trace peaks eluting on the tail of a major component.
- The entire analysis has a run time of less than 15 minutes

### Optional Configurations

- Liquid sample valves for the injection of pressurized liquid samples.
- Refinery gas analysis with trace sulfurs by FPD or SCD
- Additional boiling point column for the analysis of heavy hydrocarbons (C1–C30)
- Standard analysis with the addition of trace CO by methanizer
- Custom analyzer for performing ASTM D2163, ASTM D2712, and ISO 7941
- High temperature injection for heavy fractions
- High temperature reactor effluent with percent level water
- TCD/TCD/MSD for the analysis of reactor effluent gases



1 C7+ Backflush	11 n-Butane	21 3-Methyl-1-Butene
2 Methane	12 t-2-Butane	22 t-2-Pentene
3 Ethane	13 1-Butane	23 2-Methyl-2-Butene
4 Ethylene	14 Isobutene	24 1-Pentene
5 Propane	15 c-2-Butane	25 2-Methyl-1-Butene
6 Cyclopropane	16 Neopentane	26 c-2-Pentene
7 Propylene	17 Isopentane	27 Neohexane
8 Acetylene	18 Methyl acetylene	28 Hexane
9 Isobutane	19 Isobutane	29 Heptane
10 Propadiene	20 1,3-Butadiene	30 Benzene

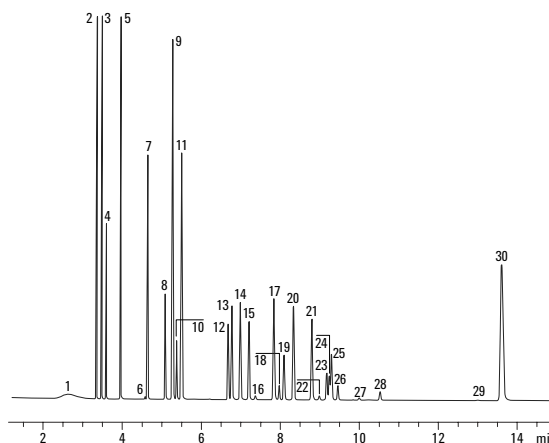

### For More Information

For more information on our products and services, visit our Web site at [www.agilent.com/chem](http://www.agilent.com/chem).




Agilent Technologies



WASSON-ECE  
INSTRUMENTATION

1 C7+ Backflush	11 n-Butane	21 3-Methyl-1-Butene
2 Methane	12 t-2-Butane	22 t-2-Pentene
3 Ethane	13 1-Butane	23 2-Methyl-2-Butene
4 Ethylene	14 Isobutene	24 1-Pentene
5 Propane	15 c-2-Butane	25 2-Methyl-1-Butene
6 Cyclopropane	16 Neopentane	26 c-2-Pentene
7 Propylene	17 Isopentane	27 Neohexane
8 Acetylene	18 Methyl acetylene	28 Hexane
9 Isobutane	19 Isobutane	29 Heptane
10 Propadiene	20 1,3-Butadiene	30 Benzene

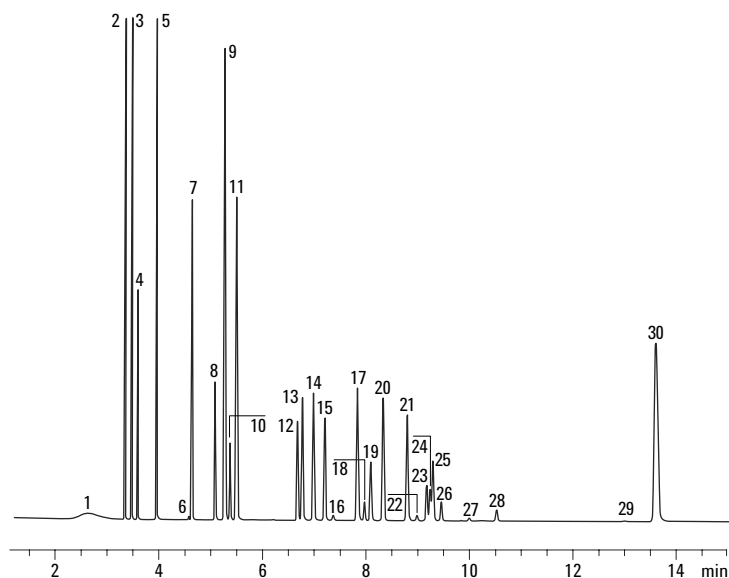

FID output from the Agilent Refinery Gas Analyzer.

Agilent shall not be liable for errors contained herein or for incidental or consequential damages in connection with the furnishing, performance, or use of this material.

Information, descriptions, and specifications in this publication are subject to change without notice.

© Agilent Technologies, Inc. 2002

Printed in the USA  
December 6, 2002  
5988-7670EN



Agilent Technologies



WASSON-ECE  
INSTRUMENTATION